



Vaporized Hydrogen Peroxide Sterilization

BY LARRY TALAPA, SENIOR SPECIALIST APPLICATION ENGINEER—3M MEDICAL SOLUTIONS DIVISION

This series of self-study lesson plans was developed by HSPA. Certified Registered Central Service Technician (CRCST) lessons provide members and certification holders with ongoing education on core Sterile Processing-related topics and processes. Purdue University's Extended Campus and HSPA both offer grading opportunities.

Earn Continuing Education Credits

Online: Visit www.myhspa.org for online grading.

By mail: For written grading of individual lessons, send completed quiz and \$15 to: Purdue Online – Central Service Young Hall, Room 405 155 S. Grant St. West Lafayette, IN 47907

Scoring: Each quiz with a passing score is worth two continuing education (CE) credits toward your CRCST recertification (12 credits).

Subscription series: Purdue Extended Campus offers an annual mail-in or online self-study lesson subscription for \$75 (six specific lessons worth two credits each toward CRCST recertification of 12 credits). Call 765.494.0959 for details.

More information: HSPA provides online grading service for any of the Lesson Plan varieties. **Purdue University provides grading services ONLY for CRCST and CIS lessons.** Direct any questions about online grading to HSPA at 312.440.0078. Questions about written grading are answered by Purdue University at 765.494.0959.

LEARNING OBJECTIVES

1. Describe the change in complexity of vaporized hydrogen peroxide sterilization in healthcare facilities since its initial use in the early 1990s
2. Recognize significant variables and practices that affect vaporized hydrogen peroxide sterilization in healthcare facilities
3. Discuss best practices for the successful use of vaporized hydrogen peroxide sterilization in healthcare facilities

The first use of vaporized hydrogen peroxide (VH₂O₂) sterilization in U.S. healthcare facilities was in 1993. At that time, the sterilizer had one cycle, one injection of VH₂O₂ sterilant, and a very limited number of compatible devices and packaging types; however, this sterilizer was a brand-new technology for the industry.

Objective 1: Describe the change in complexity of vaporized hydrogen peroxide sterilization in healthcare facilities since it was introduced in the early 1990s

Twenty-five years after its introduction, the inaugural sterilizer became obsolete and is no longer supported by the manufacturer. Today, there are multiple VH₂O₂ sterilizer manufacturers and sterilizer models, and over 20 different VH₂O₂ sterilization cycles in the U.S. market. These sterilizers use different technologies, the cycles have different sterilant injection numbers, sterilant


exposure times, VH₂O₂ concentration levels, and cycle pressure profiles. **Table 1** summarizes some of the differences between VH₂O₂ sterilization used today versus 30 years ago.

Table 1 features VH₂O₂ cycle pressure graphs and illustrates the dramatic change in the methods employed for VH₂O₂ sterilization over the last 30 years.^{2,16} A cycle pressure graph helps to illustrate the mechanism the sterilizer utilizes for sterilization. The left pressure graph in **Diagram 1** is the mechanism of the first healthcare sterilizer utilizing VH₂O₂.² As you can see from the graph, the pressure profile is very similar to the stages in today's steam and ethylene oxide (EO) sterilization cycles (e.g. air removal, sterilant injection, sterilant hold (exposure), and sterilant removal (plasma in this example)). As our understanding of VH₂O₂ sterilization developed, the VH₂O₂ cycles of today have become much more complex as depicted in the one example in the graph on the right side of **Diagram 1**.¹⁶



VH2O2 Sterilizers and Sterilization Cycles in the U.S. Healthcare Facilities ¹⁻⁹							
Year	Number of Manufacturers	Number of Sterilizer Models	Number of Sterilizer Cycles	Number of Sterilant Injections Per Cycle	Estimated VH2O2 Sterilant Concentration (mg/L)	Total Sterilant Exposure Time (min.)	Estimated Total Cycle Time (min.)
1993	1	1	1	1	6	50	> 75
2023	4*	10+	20+	2-4	6-96.6*	6-32	16-60

*Includes VH2O2 plus ozone sterilizer⁹
 Table 1

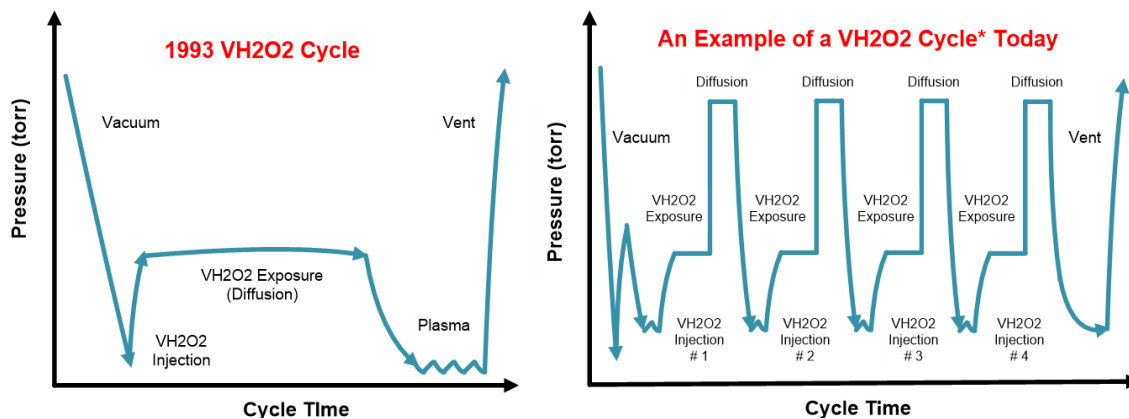

*STERIS V-PRO™ Jumen cycle¹⁶
 Diagram 1: Pressure graphs illustrating the change in VH2O2 sterilizer cycle complexity

Device compatibility for VH2O2 sterilization has grown from simple devices, like batteries and relatively wide-channeled devices, to laparoscopic instruments with narrow metal channels (0.7mm inner diameter), long, single-channel flexible endoscopes (1050 mm in length), large endoscopes for the most advanced robotic instrumentation (8.9 pounds in weight for the endoscope and sterilization tray) and multi-channel flexible endoscopes (3500 mm in length, indicated for the hydrogen peroxide plus ozone sterilizer only).⁴⁻⁹

In 2016, the U.S. Food and Drug Administration (FDA) cleared the first rapid readout biological indicator (BI) for VH2O2 sterilization. This new BI was developed with the same rapid readout technology used in BIs to

monitor steam and ethylene oxide (EO) sterilization processes for the last 20 years. A year or more later, two more rapid readout BIs for VH2O2 received FDA clearance.

All three of these new VH2O2 BIs use the same principal rapid readout technology, provide a final result in just a few minutes, and most notably, present an increased challenge to the VH2O2 sterilization process compared to the conventional readout BIs used over the last 20 years. The decreased BI readout time and increased challenge the new BIs present to the VH2O2 process, combined with the increased complexity of the VH2O2 sterilizers, cycles and load items, has led to the profession's increased awareness of the technique sensitivity for VH2O2. This awareness

sparked a culture change in the U.S. to correct and stay the course for VH2O2 sterilization.

Objective 2: Recognize significant variables and practices that effect vaporized hydrogen peroxide sterilization in health care facilities

VH2O2 sterilization is technique sensitive, which describes the variability introduced by the end user or VH2O2 sterilizer operator that can have a significant impact on the outcome of the VH2O2 sterilization process. The majority of VH2O2 sterilization processes provide a set/fixed amount of sterilant for each cycle type and load placed in the chamber. There are no make-ups of sterilant during the

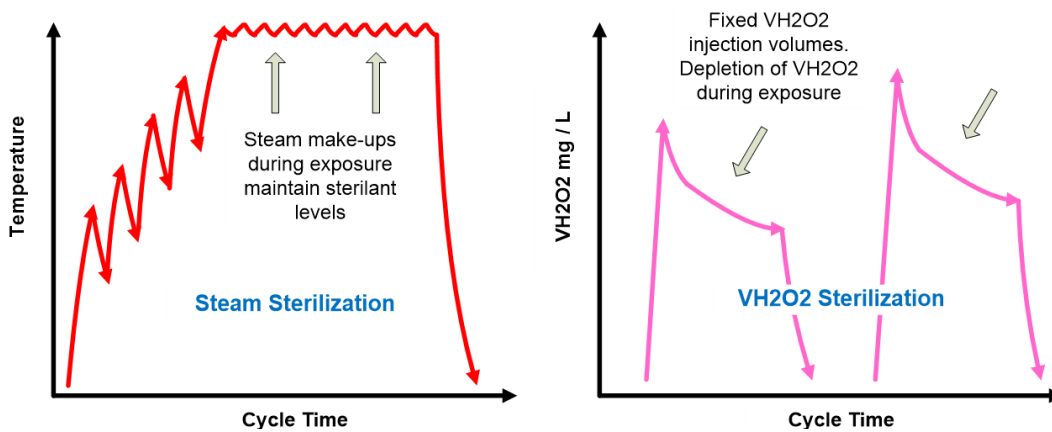


Figure 1: Relative sterilant levels in steam and VH2O2 sterilization

sterilant exposure phase; therefore, a very small or large load is exposed to the same amount of sterilant during the exposure phase. Each VH2O2 cycle can be compared to an oven that only has one temperature setting for every recipe.

The fixed amount of VH2O2 injected is relatively unstable and readily depletes during the exposure phase via several different chemical mechanisms.^{1,10,11}

Figure 1 illustrates the relative differences in sterilant levels maintained in steam sterilization versus the natural depletion that occurs during exposure after VH2O2 injection. Today, VH2O2 sterilizers are validated and cleared by the FDA, with a maximum weight limit for individual loads and each cycle type. **Table 2** shows weight limits for each VH2O2 sterilizer model and cycle.⁴⁻⁹ *Note: Exceeding the weight limit for the load can result in an automatic cycle cancellation and/or failure of quality monitoring tools. Always refer to the sterilizer manufacturer's instructions for use for specific restrictions on devices allowed for each cycle type.*

VH2O2 sterilization processes are not compatible with excessive moisture in and around devices and packaging. Excess moisture can cause automatic cycle cancellations and failure of quality

monitoring tools, resulting in rejected sterilization cycles.¹³ Temperature is a critical process parameter for VH2O2 sterilization.²⁰ This includes the temperature of the devices, packaging, and the environment of the Sterile Processing Department. The temperature of the load and department where the VH2O2 sterilizer is installed can have a negative impact on the process. If the temperature is too cool, excessive condensation of the fixed amount of VH2O2 sterilant can occur.^{6,7}

Materials compatibility is also essential for successful VH2O2 sterilization. Use of the wrong materials could result in a dramatic failure of the process. All materials placed in a VH2O2 sterilization process will affect the relatively unstable VH2O2 molecule in some way, but the user must be aware that some materials (e.g. some plastics versus some metals) can have a much more dramatic effect on the available VH2O2 by absorbing, adsorbing or decomposing VH2O2 at a higher rate.^{1, 10-12}

The use of extra (nonessential) materials in VH2O2 sterilization is another variable that is dependent on the user and can introduce significant variation to the VH2O2 sterilization

process. For example, foam tray liners, polyethylene sheet tray liners, underneath guard liners, bubble wrap tray liners and tray protectors, rubber corner protectors, foam pocketed instrument protectors, CI indicator holders, transport trays, oversized disposable sterilization wrap, 600- and 650-weight disposable sterilization wrap, and preformed disposable wraps are all examples of extraneous or nonessential materials used in healthcare facilities. Again, because VH2O2 cycles use a fixed amount of sterilant, best practice would be to limit or eliminate the use of any extra materials that could absorb the fixed amount of available VH2O2 sterilant.

ANSI/AAMI ST58:2013 (R2018) *Chemical sterilization and high-level disinfection in health care facilities*¹⁷ and AORN Guidelines for Perioperative Practice¹⁸ are our standard references in the U.S. for the use of VH2O2 sterilization in healthcare facilities. Both references point to some of the items discussed previously in this lesson but are not explicit on many items that help ensure a successful VH2O2 sterilization cycle. ANSI/AAMI ST58:2013 is currently under revision by AAMI Working Group 61.



Chamber weight limits per common sterilizer model and cycle types			
	Model	Cycle	Weight (lb.) Limit
Advanced Sterilization Products (ASP) [®]	STERRAD [®] 100S	Standard (default) ⁵	Not defined
	STERRAD [®] NX	STANDARD ¹⁴	10.7
		ADVANCED ¹⁴	10.7
	STERRAD [®] 100NX	STANDARD ¹⁵	21.4
		FLEX ¹⁵	21.4
		EXPRESS ¹⁵	10.7
		DUO ¹⁵	13.2
STERIS [®]	V-PRO maX 2	Non-Lumen ⁸	50.0
		Lumen ⁸	19.6
		Flexible ⁸	24.0
		Fast Non-Lumen ⁸	11.0
STYRKER [®]	STERIZONE [®] VP4* (VH2O2 plus ozone sterilizer)	Cycle ¹⁹	75.0

Table 2

Objective 3: Discuss best practices for the successful use of vaporized hydrogen peroxide sterilization in healthcare facilities

Following the device manufacturer’s instructions for use (IFU) seems straight forward and unassuming, but one may be surprised by what can be uncovered verifying each detail in the IFU for each device the facility sterilizes using VH2O2. Let’s explore a common scenario regarding failed cycles.

For the subject device, Intuitive Surgical’s da Vinci Xi[®] endoscope processed in the ASP[®] STERRAD[®] 100NX EXPRESS cycle, the length of the endoscope is approximately 600 mm, and the diameter of the shaft is 8.8 mm. The maximum weight of the tray and endoscope is 8.9 pounds. This represents of the largest devices in the U.S. labeled for VH2O2 sterilization. The ASP[®] STERRAD[®] 100NX EXPRESS cycle is validated for a load with a maximum

weight limit of 10.7 pounds, loaded only on the sterilizer chamber’s bottom shelf. This cycle has the shortest total VH2O2 exposure time (estimated six minutes) for any VH2O2 sterilization cycle currently on the U.S. market. When combining of the largest devices with the shortest total VH2O2 sterilitant exposure time, it is imperative to diligently follow the IFU for the da Vinci Xi endoscope and ASP[®] STERRAD[®] 100NX EXPRESS cycle to ensure consistent, successful cycles. The IFU for the da Vinci Xi endoscope processed in the STERRAD[®] 100NX EXPRESS cycle notes the following¹⁹:

- Confirm endoscope is properly loaded into tray (PN 400490)
- Do not stack trays during sterilization
- Do not process more than one tray at a time
- Only process one tray on the bottom shelf
- Only use the Express Cycle (on the STERRAD[®] 100NX)

- Always use sterilization wrap rated for:
 - » 9-13 pounds (medium weight or lighter) or
 - » 400-weight thickness or lighter
- for the plastic tray (PN 400490) with the STERRAD 100NX (Express Cycle)
- Using a thicker wrap may result in incomplete sterilization of the Xi Endoscope

Again, moisture is not compatible with VH2O2 sterilization. Ensuring the da Vinci Xi endoscope is dry and verifying there is no water trapped in the device (button flush ports, input discs, housing and plastic tray) is crucial to ensuring successful VH2O2 sterilization. Verifying all recommended drying steps are completed per the IFU can solve and prevent failed VH2O2 sterilization cycles in the previously described scenario.

This saying rings true for a best-practice VH2O2 sterilization: know what you’re loading; load only what you know. It is imperative that all operators of VH2O2 sterilizers understand the composition of each load placed in the sterilization chamber. Some basic questions for this best practice include:

- Are the devices labeled for their specific VH2O2 sterilizer model and cycle type?
- Is the total load weight below the validated and cleared weight limit?
- Is the packaging type acceptable for use in VH2O2, and is the device weight under the limit for the packaging type?
- Could the device be labeled for another sterilization method (such as steam)?
- Are there any non-essential extraneous packaging items that could be avoided?
- What is the total material composition of the load? Is the load overly weighted with items that have a higher propensity of depleting the fixed amount of VH2O2?



Understanding these basic variables for each VH₂O₂ load will help the user better understand the effect these factors have on the process and will help ensure consistent, successful process outcomes.

Many SP technicians have worked in more than one facility in the same city or even in several facilities across the country, which means they have experienced different processes and practices. There are correlations to failed VH₂O₂ sterilization cycles based on practices brought from one facility to another facility by a new SP technician. Providing thorough training and competency evaluations for every new employee, based on the hiring facility's policies and procedures, helps promote process and practice consistency and better outcomes. Each facility has common sets of devices routinely sterilized using VH₂O₂ as well as procedures and methods for drying devices, packaging devices, loading devices. Further, they may use non-essential (extraneous) packaging items and different sterilizer models and cycle types. Because of the technique sensitivity of VH₂O₂ sterilization, a small shift in procedures by a new SP technician may result in a sporadic failed VH₂O₂ sterilization cycle.

Another commonly observed occurrence is the increased frequency of BI monitoring for every load, combined with quarantining the load until the BI result is known. ANSI/AAMI ST58, Section 9.5.4.3, states that "...[a] BI should also be used at least daily, but preferably in every sterilization cycle."¹⁷ The Association of periOperative Registered Nurses Guideline for Sterilization, Section 10.10.2) is slightly more specific: "Perform routine sterilizer efficacy monitoring every day the sterilizer is used for each cycle type ... preferably with each load."¹⁸ In hospitals, end users typically place

a BI and internal CI in a peel pouch designed for use in VH₂O₂ sterilizers and position the pouched BI in the sterilizer chamber as recommended by the sterilizer manufacturer. When users switch to a new BI that provides a result in minutes versus days, they can quickly move to every load monitoring to provide a consistent level of patient care. In addition, the same users now quarantine every VH₂O₂ load until the BI result is known to mitigate the risk of large recalls in the event of a cycle failure.

Unfortunately, misinformation exists regarding the use of BIs for VH₂O₂ sterilization. Because an international standard does not yet exist, the global healthcare industry has no standardization on performance requirements for BIs used in VH₂O₂. In the U.S., BIs are medical devices; the FDA regulates BIs used in healthcare facilities and has a set of testing requirements for VH₂O₂ BIs cleared for use in the U.S. The FDA is the highest authority in the U.S. (not the sterilizer manufacturer) on the final decision on which BIs are cleared as compatible (safe and effective) for use in healthcare facilities' VH₂O₂ sterilizers. Many users were unaware that there is no requirement for a sterilizer manufacturer to validate nor endorse indicators designed to monitor their sterilizers. The decision regarding the safety and efficacy of sterilization monitors is addressed by the FDA's review and clearance procedures. There are many examples of monitoring products from multiple manufacturers being used to monitor steam, EO and VH₂O₂ sterilizers.

Conclusion

Many factors contribute to safe, effective and efficient VH₂O₂ sterilization, including user knowledge about their

facility's policies, procedures and practices and consistent adherence to the manufacturers' IFU. Ensuring safe, effective VH₂O₂ outcomes also requires thorough training for SP technicians and strict adherence to industry standards, guidelines and best practices.

REFERENCES

1. Peroxygen Compounds. *Block's Disinfection, Sterilization, and Preservation*, fifth ed., Chapter 9. 2001.
2. Russell, Hugo and Ayliffe's *Principles and Practice of Disinfection, Preservation and Sterilization*, fourth ed., Chapter 12. 2004.
3. STERIS® Document# M3644EN2012-10, Rev. B
4. FDA 510(k) K991999 _ STERRAD® 100S STERILIZER
5. FDA 510(k) K111377 _ STERRAD® 100NX STERILIZER DUO CYCLE
6. FDA 510(k) K042116 _ STERRAD® NX Sterilizer
7. FDA 510(k) K071385_ STERRAD® 100NX Sterilizer
8. FDA 510(k) K172754_V-PRO® maX 2 Low Temperature Sterilization System
9. STERIZONE® VP4® Low Temperature Sterilizer Brochure 04/2018
10. Kohler et al. U.S Patent 6,528,016 B1. March 4, 2003.
11. McVey U.S Patent 6,875,399 B2. April 5, 2005.
12. Childers et al. U.S Patent 5,830,409. November 3, 1998.
13. Robinson N, Eveland R. "Using HPG sterilization for heat-sensitive devices." *Healthcare Purchasing News*, January 2015.
14. FDA 510(k) K160818 _ STERRAD® NX Sterilizer with ALLClear™ Technology
15. FDA 510(k) K160903_ STERRAD® 100NX Sterilizer with ALLClear™ Technology
16. Technical Data Monograph STERIS V-PRO maX Document #10037127 Rev. C (M3635EN)
17. ANSI/AAMI ST58:2013 (R2018) *Chemical sterilization and high-level disinfection in health care facilities*
18. AORN *Guidelines for Perioperative Practice*. 2023.
19. Intuitive_Surgical_da_Vinci_Surgical_System_Prep_Pack_554063_02_rB_2019
20. ISO 11140-1:2014 *Sterilization of healthcare products – Chemical indicators – Part 1: general requirements*



CRCST Self-Study Lesson Plan Quiz

Vaporized Hydrogen Peroxide Sterilization

Lesson No. CRCST 190 (Technical Continuing Education – TCE) · Lesson expires May 2026

- Vaporized hydrogen peroxide (VH2O2) sterilization was first introduced to healthcare facilities:
 - With limited cycles and packaging selections
 - With multiple cycles and packaging selections
 - With multiple injections of hydrogen peroxide
 - In the early 2000s
- A cycle pressure graph illustrates:
 - The process flow for all sterilization methods
 - Device compatibility for VH2O2 cycles
 - The mechanism the sterilizer uses for sterilization
 - When the sterilizer is due for maintenance
- The pressure profile of VH2O2 sterilizers is:
 - Similar to the first hydrogen peroxide sterilizers
 - Similar to ethylene oxide and steam sterilizer cycles
 - Representative of the total chamber weight
 - None of the above
- Today's VH2O2 sterilizers can sterilize:
 - All items requiring low-temperature sterilization
 - Some multi-channeled flexible endoscopes, laparoscopic devices and robotic devices
 - Only wide-channeled devices
 - None of the above
- VH2O2 is technique sensitive because of the:
 - Variability introduced by the sterilizer operator
 - Varying temperatures during the cycle
 - Lack of instructions for use
 - All the above
- When using VH2O2, it is important to monitor:
 - Moisture on the instruments
 - The materials being processed in the sterilizer
 - Temperature of the devices and the environment
 - All of the above
- As long as the chamber is loaded properly, exceeding the validated load maximum weight limit has no effect on the sterilization process.
 - True
 - False
- Which of the following materials can affect the VH2O2 sterilization process?
 - Rubber corner protectors
 - Oversized sterilization wrappers
 - Foam instrument protectors
 - All the above
- It is considered best practice for VH2O2 sterilization to follow:
 - The Centers for Disease Control and Prevention's low-temperature guidance document
 - ANSI/AAMI ST58 and Association of perioperative Registered Nurses *Guidelines for Perioperative Practice*
 - The U.S. Food and Drug Administration's (FDA's) Good Manufacturing Practice document
 - Safety protocols set forth by the Occupational Safety and Health Administration
- Which of the following is a best practice for loading VH2O2 sterilizers?
 - Knowing what to load
 - Not exceeding the weight load for the sterilizer
 - Using the correct package for the sterilizer
 - All the above
- Sterile Processing technicians with experience using VH2O2 sterilization at another facility still need competency assessments performed at the new facility.
 - True
 - False
- VH2O2 cycles should be monitored with a biological indicator (BI):
 - Monthly or, preferably, weekly
 - Weekly or, preferably, daily
 - At least daily
 - At least daily but preferably for every cycle
- Many facilities are increasing the frequency of:
 - Monitoring for hydrogen peroxide concentrations
 - Monthly competency evaluations for all staff members
 - Biological monitoring for VH2O2 sterilization cycles
 - VH2O2 chamber cleaning
- Which of the following is considered the highest authority for the clearance of BIs in the US?
 - Occupational Safety and Health Administration
 - FDA
 - Sterilizer manufacturer
 - BI manufacturer
- Users must approach the use of VH2O2 with:
 - Great discipline
 - Knowledge of the FDA's labeling document
 - Knowledge of the sterilizer's manufacturing process
 - All the above

Request for paper/pencil scoring

(Please print or type information)

- I have enclosed the scoring fee of \$15. (Please make checks payable to Purdue University. We regret that no refunds can be given.)
- Check here if you have a change of address.
- Check here if you wish to have your results emailed to you.

Detach quiz, fold, and return to:

Purdue Online – Central Service
 Young Hall Room 405
 155 S. Grant St.
 West Lafayette, IN 47907

Name _____

Mailing Address (Be sure to include apartment number or post office box.) _____

City _____

State/Province _____

Zip Code/Postal Code _____

Daytime Telephone _____

HSPA Membership Number _____

Email Address _____

If your name has changed in the last 12 months, please provide your former name. Purdue University is an equal access/equal opportunity institution.